


# VariLoc

## Solutions to Foot & Ankle Fractures

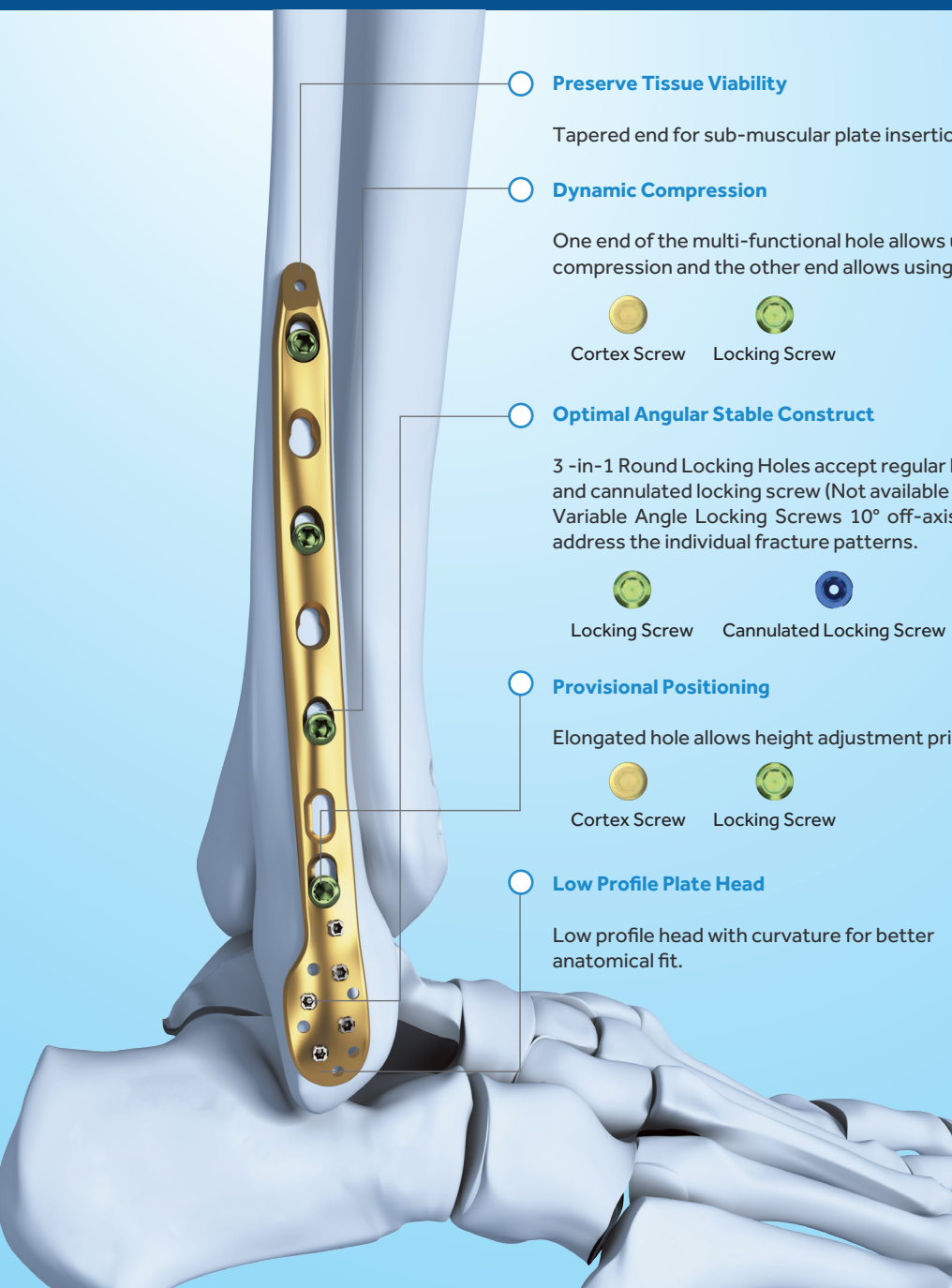

### Preserve Tissue Viability

Tapered end for sub-muscular plate insertion which may help to preserve tissue viability.

### Dynamic Compression

One end of the multi-functional hole allows using cortex screw to perform dynamic compression and the other end allows using locking screws for angular stability.



Cortex Screw



Locking Screw

### Optimal Angular Stable Construct

3 -in-1 Round Locking Holes accept regular locking screw, Variable Angle (VA) locking screw and cannulated locking screw (Not available for 2.7mm system).

Variable Angle Locking Screws 10° off-axis screw angulation in all directions in order to address the individual fracture patterns.



Locking Screw



Cannulated Locking Screw



Variable Angle Locking Screw

### Provisional Positioning

Elongated hole allows height adjustment prior to the insertion of the other screws in the shaft.




Cortex Screw



Locking Screw

### Low Profile Plate Head

Low profile head with curvature for better anatomical fit.

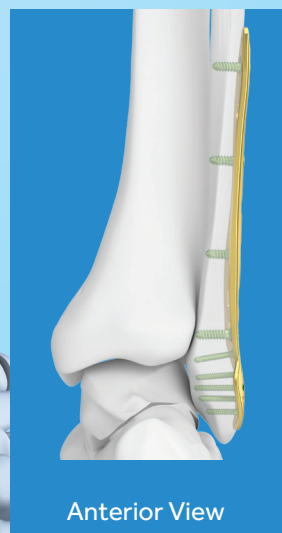

Anterior View



Inferior View

# PLATE AND SCREW OPTIONS

Description	Image	Plate Size	Length (mm)	Screw	
Lateral Distal Fibular Locking Compression Plate 2.7/3.5		5H Head 3-15H Shaft	86-242	Head	2.7 locking screw 2.7 VA locking screw
				Shaft	3.5 cortex screw 3.5 locking screw
Anterolateral Distal Tibia L Locking Compression Plates, 3.5		4H Head 5-15H Shaft	117-277	Head	3.5 VA locking screw 3.5 locking screw 4.0 cannulated locking screw
				Shaft	3.5 locking screw 3.5 cortex screw
Lateral Distal Fibular Locking Compression Plate II 2.7/3.5		5H Head 3-15H Shaft	86-242	Head	2.7 locking screw 2.7 VA locking screw
				Shaft	3.5 cortex screw 3.5 locking screw
Distal Tibia Locking Compression Plates II, 2.7/3.5		8H Head 4-14H Shaft	116-246	Head	3.5 locking screw 3.5 VA locking screw 4.0 cannulated locking screw
				Shaft	3.5 locking screw 3.5 cortex screw
Distal Medial Tibia Locking Compression Plates (with tip), 2.7/3.5		9H Head 4-14H Shaft	116-246	Head	2.7 VA locking screw 2.7 locking screw
				Shaft	3.5 locking screw 3.5 cortex screw
Calcaneus Locking Plate II 3.5		13-16H	60-78	3.5 locking screw	
Calcaneus Locking Plate III 3.5		12H	60-68	3.5 locking screw	